TUNTURI®

M40 MASSAGE GUN

GB	Assembly manual	07 - 11
DE	Montageanleitung	12 - 16
FR	Manuel d'assemblage	17 - 21
NL	Montagehandleiding	22 - 26


Attention, Achtung, Attention, Attentie

- **Please read this assembly manual carefully prior to assembling this product.**
- Lesen Sie sich vor der Montage des Produkts diese Montageanleitung sorgfältig durch.
- Veuillez lire attentivement le présent manuel d'assemblage avant de monter ce produit.
- Lees deze montagehandleiding zorgvuldig voordat u dit product in elkaar gaat zetten.


A

B

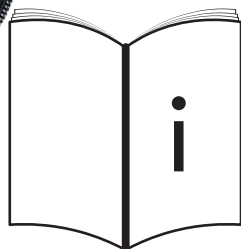C

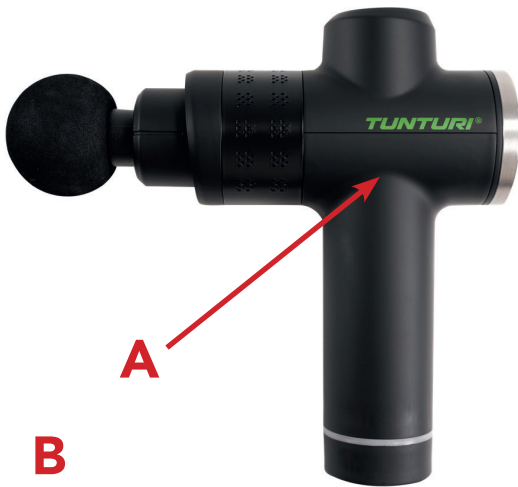B

C

D

E

F

G

M40 Massage Gun

D-1

D-2

D-3

D-4

D-5

D-6

USER MANUAL ON OUR WEBSITE

To reduce the impact on the environment, you will only find the assembly steps in this manual. To learn how to use and maintain this product, please download the user manual from:

manuals.tunturi.com

Index

Safety warnings.....	7
Usage instructions	8
Charging.....	8
Operating	8
Button functions.....	8
Maintenance, cleaning & storage	9
Troubleshooting	9
Error codes.....	9
Defects & malfunctions	9
Specifications.....	10
Massage head user instructions	10
Declaration of the manufacturer	11
Disclaimer	11

Massage Gun

Safety warnings

⚠ WARNING

- Read the safety warnings and the instructions. Failure to follow the safety warnings and the instructions can cause personal injury or damage to the equipment. Keep the safety warnings and the instructions for future reference.

⚠ DANGER

- To reduce risks of electric shock, fire and personal injury, this product must be used in accordance with the following instructions.
- The equipment is suitable for domestic use only. The equipment is not suitable for commercial use.
- The use of this equipment by children or people with a physical, sensory, mental or motorial disability, or lack of experience and knowledge can give cause to hazards. People responsible for their safety must give explicit instructions or supervise the use of the equipment.
- If you experience nausea, dizziness or other abnormal symptoms, immediately stop your workout and consult a physician.
- The equipment is suitable for indoor use only. The equipment is not suitable for outdoor use.
- Only use the equipment in environments with ambient temperatures between 10°C and 35°C. Only store the equipment in environments with ambient temperatures between 5°C and 45°C.
- Do not use or store the equipment in humid surroundings. The air humidity must never be more than 80%.

- Only use the equipment for its intended purpose. Do not use the equipment for other purposes than described in the manual.
- Do not use the equipment if any part is damaged or defective. If a part is damaged or defective, contact your dealer.
- Keep your hair away from the moving parts.
- Keep jewellery and other objects away from the moving parts.
- Keep fingers, hair or any other body part away from the shaft and backside of the massage head as pinching may occur.
- Make sure that only one person uses the equipment at a time.
- Do not open equipment without consulting your dealer.
- Use only on a dry clean surface of the body by lightly pressing and moving across the body for approximately 60 seconds per region.
- Use the Massage Gun only on the body's soft tissue as desired without producing pain or discomfort.
- Do not use on head or any hard or bony area of the body.
- Use only the massage heads that provide the best desired effect.
- Bruising can occur regardless of control setting or pressure applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort
- Do not place any objects into the ventilation holes of the Massage Gun.
- Do not immerse in water or allow water to enter the ventilation holes of the Massage Gun.
- Do not drop or misuse the Massage Gun.
- Only recharge the massage gun with a 5 volt USB charger in combination with an USB-C cable.
- Carefully examine device and battery before each use.
- Do not tamper or alter Massage Gun in any way.
- Never leave the Massage Gun operating or charging unattended.
- Keep ventilation ports free from dust and debris.
- Do not operate continuously for more than one hour. Allow device to rest for 30 minutes before re-using.
- Unplug the charger from device battery after charging or prior to use.
- Do not use if injured. Consult your doctor before using this product.

Usage instructions

Charging

- Fully charge battery for at least 3 hours before first use.
- To charge, connect the USB cable type-C end of the supplied cable to the USB charging port on the massage gun handgrip's bottom end. Connect the USB cable Type-A end of the supplied cable to a USB charger and plug the charger into a power socket.
- LED ring on handgrip lights pulsate to show battery level and indicate active charging.
- Full charge is indicated when the LED ring light stops flashing and remains fully illuminated.
- Battery may be recharged at any time and at any battery level.
- It is not recommended to fully deplete the battery to the red LED level.

!! NOTE

- When the battery power is less than 20%, the light ring at the bottom of the machine shows lights red.

⚠ WARNING

- We strongly recommend to not operate the machine during charging.

Operating

- Power device on by moving the power switch on the bottom of the battery to the ON position when being in OFF mode (Fig. C - C)
- Massage desired body part at the desired speed applying pressure as needed.

Button functions

Power Switch (Fig. C - C)

Turn On/ Off the power of the machine by setting the switch in desired position.

On/ Off button (Fig. C - F)

- In Off mode, press button once to start massage function.
- In On mode, press button once to stop massage function.

+/- Button (Fig. C - E&G)

Press + and - buttons to adjust machine speed.

!! NOTE

- If button not respond, you have reached your minimum/ maximum speed.

Maintenance, cleaning & storage

- Use a slightly damp towel to clean device and a dry with a soft cloth.
- Avoid moist to enter the device's housing at all times.
- For storage or travel store in supplied storage box.

!! NOTE

- When not using the massage gun for a longer period of time, we recommend you to charge the device at least once a month. This is to keep the battery in optimal condition.
- An exhausted battery will shorten it's lifetime or cause performance failure.

Troubleshooting

When there is too much stress applied to the massage gun when being in use, the massage gun will shut off.

Try to resolve by:

- When the massage gun is hot, give it some time to cool off.
- Press the control button and hold for 3 seconds to turn in Stand By mode again, and operate it as normal.

When a abnormal sounds comes out of the massage gun during usage, stop using the massage gun to avoid damage.

Try to solve by:

- Reinstalling the massage head to make sure the fixture is correct.
- Applying less load on the massage gun when using.
- When the massage gun will not turn on, try to charge it for at least 4 hours to make sure the battery is charged.

Error codes**E1 Battery low**

Problem: Battery power is too low to operate.

Solve by: Charge battery for a full cycle at once.

E3 Over pressure protection

Problem: Too much pressure on the motor during operation.

Solve by: Cycle Off/ On to rest the error code.

Use less pressure during operation to avoid the over pressure protection to appear again.

Defects & malfunctions

Despite continuous quality control, defects and malfunctions caused by individual parts can occur in the equipment. In most cases, it will be sufficient to replace the defective part.

- If the equipment does not function properly, immediately contact the dealer.
- Provide the model number and the serial number of the equipment to the dealer. State the nature of the problem, conditions of use and purchase date.

Specifications

The Tunturi Massage Gun is a portable rechargeable massage device with large torque and long adjustment. It can effectively reduce the pain caused by the accumulation of lactic acid after intense exercise, and relieve the discomfort of muscles in the sedentary and long-lying population.

Model	Massage Gun MG40
Dimensions (L x W x H)	15.1 x 5.6 x 18.9 cm 5.9 x 2.2 x 7.4 inch
Product weight	0.6 kg / 1.3 lb
Battery capacity	1500 mAh
Power	15W
Charger power output	5V2A
Amplitude	8mm
Stall force	14 kg
RPM	2800 rpm
Charge time	3 hours
Full charge battery use time at	
Level 20 speed	2 hours
Level 10 speed	2 hours, 50 minutes
Level 1 speed	3.5 hours

Massage head user instructions

With unit switched in OFF mode, fully insert massage head by pressing in firmly.
With unit switched in OFF mode, remove the massage head by pulling it out. To ease the release you can rotate slightly and pull at the same time.

TIP

- When remove the massage head it is best to pull and turn the massage head, as if it is a threaded connection, to make it come loose more easily.
- When place the massage head it is best to push and turn the massage head, as if it is a threaded connection, to make secure the fixation more easily.

WARNING

- Do not press/ pull the massage head in/ out at an angle to avoid damages.

Standard ball (Fig. D-1)

For core,glutes and legs.

Flat head (Fig. D-2)

For chest, abdomen and arms.

U shape (Fig. D-3)

For shoulders and spine

Conical shape (Fig. D-4)

For soles of feet and palms.

Mushroom head (Fig. D-5)

Same area as standard ball and flat head.

Wedge (Fig. D-6)

For shoulder blades and calves

Tunturi purchaser's warranty

Warranty terms

The consumer is entitled to the applicable legal rights stated in the national legislation concerning the commerce of consumer goods. This warranty does not restrict these rights. The Purchaser's Warranty is only valid if the item is used in an environment approved by Tunturi New Fitness BV for that particular equipment. The product-specific approved environment is stated in the Owner's Manual provided with your equipment.

Warranty conditions

The warranty conditions commence from the date of purchase. Warranty conditions may vary per country, so please consult your local dealer for the warranty conditions.

Warranty coverage

Under no circumstances shall Tunturi New Fitness BV or the Tunturi Distributor be liable by virtue of this warranty or otherwise for any special, indirect, secondary or consequential damages of any nature arising out of any use or inability to use this equipment.

This warranty excludes.

- Damage caused by accident, abuse, mishandling, or transport;
- Units not used in accordance with instructions;
- Damage exceeding the cost of the product;

- Deterioration of the delivered product resulting from abnormal storage and / or safeguarding conditions on the client's premises;
- Failure to provide the dated proof of purchase.

Failure to follow the instructions given in the Owner's Manual will invalidate the product's warranty.

Declaration of the manufacturer

Tunturi New Fitness BV declares that the M40 Massage Gun Massage device is in conformity with the following standards and directives:

- Low Voltage Directive (2014/35/EU)
- Electromagnetic Compatibility Directive (2014/30/EU)
- RoHS Directive (2011/65/EU)
- Battery Regulation (EU) 2023/1542.

Hereby, Tunturi New Fitness BV declares that the radio equipment type M40 Massage Gun Massage device is in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address: <https://doc.tunturi.com>

Tunturi New Fitness BV
Purmerweg 1
1311 XE Almere
The Netherlands

Disclaimer

© 2021 Tunturi New Fitness BV
All rights reserved.

- The product and the manual are subject to change.
- Specifications can be changed without further notice.
- Check our website for the latest user manual version.

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TUNTURI NEW FITNESS B.V.